


Rotary

R O T A T A B L E I B D

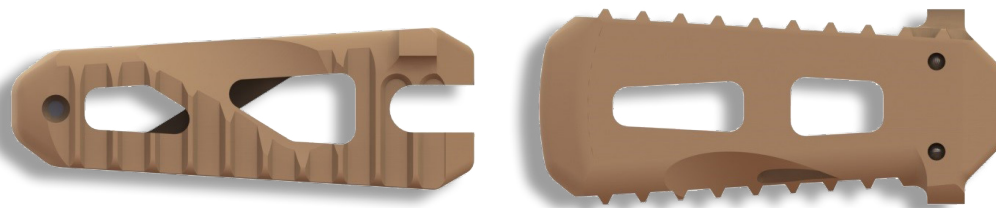

PRODUCT FEATURES:

- 🌊 Designed specifically for insert and rotate technique
- 🌊 Proprietary dual taper design aids in initial distraction and creates lordosis after rotation
- 🌊 Corner radii provide smooth 'cam' effect during rotation and distraction maneuver
- 🌊 Counter rotation barbs provide implant stability
- 🌊 Implant material: PEEK
- 🌊 Height: 9-17mm, Length: 25mm, Lordosis: 7°

ORDERING INFORMATION:

CATALOG NUMBER	PRODUCT DESCRIPTION
STL-009-006	9x25mm Lordotic Spacer, 7°
STL-011-008	11x25mm Lordotic Spacer, 7°
STL-013-010	13x25mm Lordotic Spacer, 7°
STL-015-012	15x25mm Lordotic Spacer, 7°
STL-017-014	17x25mm Lordotic Spacer, 7°